

AbsorbENT™

Nasal and Sinus Packing

RHINOLOGY



- Soft, absorbent, lint and fiber free
- Rounded edges for quick and easy insertion
- Even expansion for desired gentle mucosal pressure
- Expand up to six times their compressed size to fill cavity
- Rapid wicking of blood and fluid
- Sterile, single use, individually packaged

Our AbsorbENT™ nasal and sinus packs are three-dimensional, interconnected open cells that are extremely soft and pliable. This interconnected structure produces packing that is exceptionally strong, highly absorbent and virtually lint and fiber free.

All dimensions are in centimeters.



Epistaxis Pack with String

Product No.	Description	Unit / Box <i>Sterile, single use</i>
RH-2310-10	10.0L x 1.5W x 2.5H	10
RH-2315-20	5.5L x 1.5W x 2.5H	20



Nasal Packing

Product No.	Description	Unit / Box <i>Sterile, single use</i>
RH-7401-10	4.5L x 1.5W x 2.0H	10
RH-7404-10	8.0L x 1.5W x 2.0H	10



Nasal Packing with Airway and String

Product No.	Description	Unit / Box <i>Sterile, single use</i>
RH-7405-10	8.0L x 1.5W x 2.0H	10
RH-7408-10	4.5L x 1.5W x 2.0H	10



Nasal Packing, Anatomical with String

Product No.	Description	Unit / Box <i>Sterile, single use</i>
RH-7409-10	8.0L x 1.5W x 3.0H	10



Nasal Packing, Anatomical with Airway and String

Product No.	Description	Unit / Box <i>Sterile, single use</i>
RH-7410-10	8.0L x 1.5W x 3.0H	10



Nasal Packing, Thin

Product No.	Description	Unit / Box <i>Sterile, single use</i>
RH-7406-10	8.0L x 1.0W x 2.0H	10



Nasal Packing, Thin Anatomical

Product No.	Description	Unit / Box <i>Sterile, single use</i>
RH-7411-10	8.0L x 1.0W x 3.0H	10



Sinus Packing with String

Product No.	Description	Unit / Box <i>Sterile, single use</i>
RH-7530-10	3.5L x 0.6W x 1.2H	10*
RH-7532-10	3.5L x 0.9W x 1.2H	10*

*10 packs of 2 each per box.

Call **1.888.229.2875** to order today!

MORE



AbsorbENT™


Nasal and Sinus Packing

Doyle Patient-Friendly Nasal Packs


Developed in cooperation with Dr. Donald Doyle, MD, FACS

- Anatomically shaped, compressed sponges
- Sterile, single use, individually packaged
- Pre-cut from back to front except for web at distal end


Doyle Patient-Friendly Nasal Packs are inserted and function like a standard anatomical nasal pack. On removal, the pack is pulled slightly out of the nasal cavity to expose the uncut web. The web is then cut through with scissors, which separates the back into two parts. The two parts are removed one at a time, greatly reducing patient discomfort.

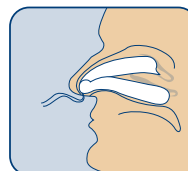

Inserted into nasal fossa.

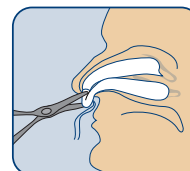

Pull out slightly and cut through remaining web.

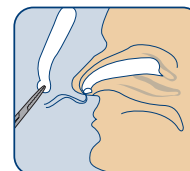

Remove two parts, one at a time.

All dimensions are in centimeters.



Anatomical Nasal Packing with String

Product No.	Description	Unit / Box <i>Sterile, single use</i>
RH-7412-10	8.0L x 1.5W x 3.0H	10



Anatomical Nasal Packing with String, Thin

Product No.	Description	Unit / Box <i>Sterile, single use</i>
RH-7414-10	8.0L x 1.0W x 3.0H	10



Anatomical Nasal Packing with Airway Tube

Product No.	Description	Unit / Box <i>Sterile, single use</i>
RH-7413-10	8.0L x 1.5W x 3.0H	10

Call **1.888.229.2875** to order today!



Summit Medical, Inc. | St. Paul, MN USA | local p: 651.789.3939 f: 651.789.3979 | toll free p: 1.888.229.2875 f: 1.888.229.1941

www.summitmedicalusa.com

CAUTION: Federal (USA) law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labeling supplied with each device.
©2011 Summit Medical, Inc. All rights reserved. ENT-0006/B APR11